

# Entrance and Exit Angle Variations

Rosalind Liang

August 2021

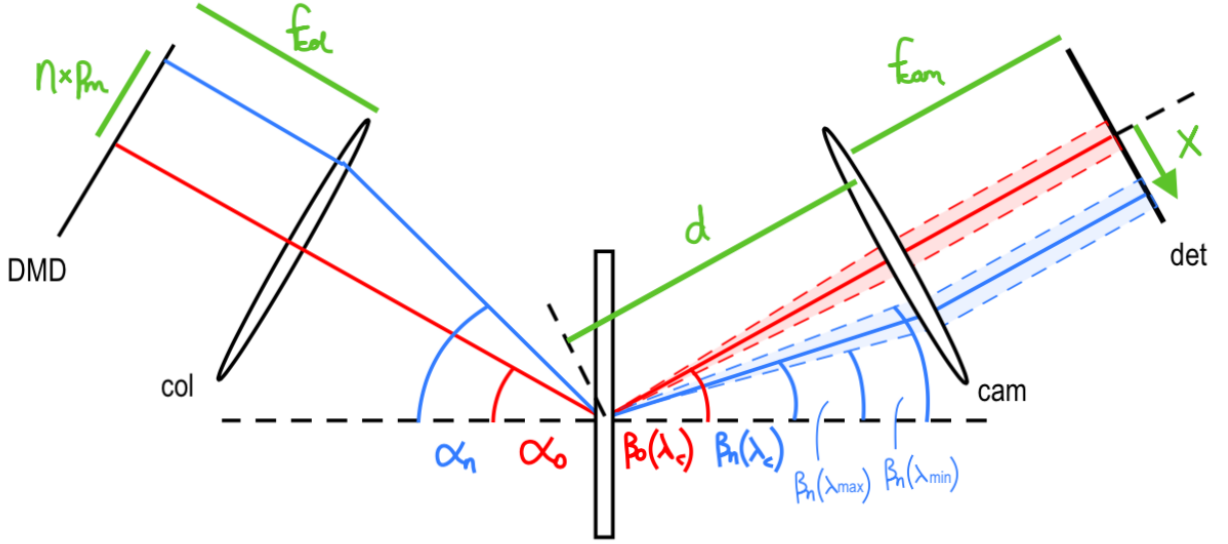


Figure 1: Schematic illustrating most of the significant variables in this analysis

Under the Bragg condition for the central wavelength  $\lambda_c$ , we can assume that the chief ray from the central slit ( $n=0$ ) has an entrance angle  $\alpha_0$  equal to its exit angle  $\beta_0(\lambda_c)$ . This angle can be calculated using the grating equation, as shown in equation 1.

$$G\lambda_c = \sin \alpha_0 + \sin \beta_0(\lambda_c) = 2 \sin \alpha_0$$

$$\alpha_0 = \arcsin\left(\frac{G\lambda_c}{2}\right) \quad (1)$$

$$= \arcsin\left(\frac{600}{mm} \times \frac{.00055mm}{2}\right)$$

$$\simeq 0.1658rad \simeq 9.497deg$$

With some trigonometry referencing figure 1, the entrance angle for slits  $n$  micromirrors from the center can be calculated using equation 2.

$$\alpha_n = \alpha_0 + \arctan\left(\frac{np_M}{f_{col}}\right) \quad (2)$$

Using the grating equation, the corresponding exit angle  $\beta_n$  can be written as a function of wavelength  $\lambda$ , as shown in equation 3.

$$\begin{aligned} G\lambda &= \sin \alpha_n + \sin \beta_n(\lambda) \\ \beta_n(\lambda) &= \arcsin(G\lambda - \sin(\alpha_n)) \end{aligned} \quad (3)$$

To calculate entrance and exit angles for the slit on the outermost edge of the DMD, we need to set  $n = n_{edge}$ . For a 2x2 micromirror slit, the value of  $n_{edge}$  is defined by equation 4 below. See figure 2 for a visualization of why this is the case.

$$\begin{aligned} n_{edge} &= \frac{N_{MX}}{2} - 1 \\ &= \frac{750}{2} - 1 = 374 \end{aligned} \quad (4)$$

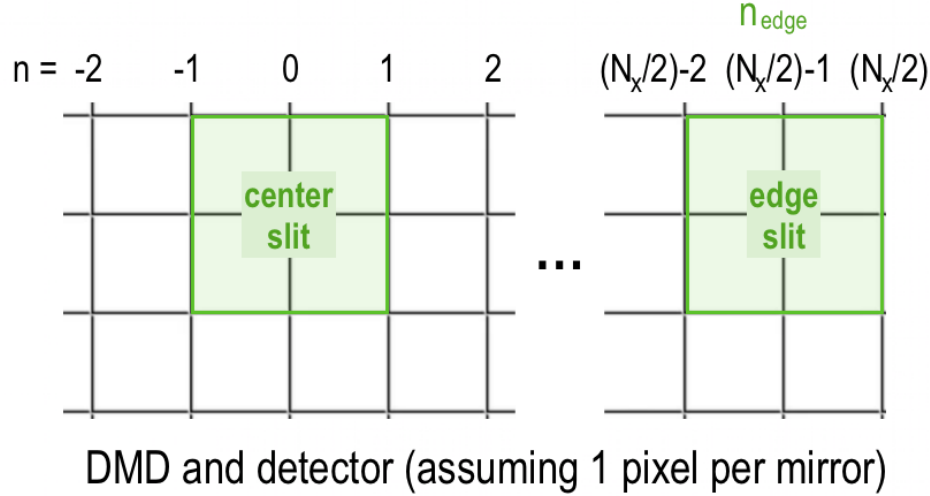


Figure 2: Diagram illustrating the allocation of slits on the micromirror/detector

Now, using equations 2 and 3, we can calculate the entrance and exit angles of the edge slit.

$$\begin{aligned} \alpha_{n_{edge}} &= \alpha_0 + \arctan\left(\frac{n_{edge}p_M}{f_{col}}\right) = 0.1658rad + \arctan\left(\frac{374(.0137mm)}{85mm}\right) \\ &\simeq 0.2260rad \simeq 12.95deg \end{aligned}$$

$$\begin{aligned} \beta_{n_{edge}}(\lambda_c) &= \arcsin(G\lambda_c - \sin(\alpha_{n_{edge}})) = \arcsin\left(\frac{600}{mm}.00055mm - \sin(.2260)\right) \\ &\simeq 0.1062rad \simeq 6.082deg \end{aligned}$$

The next step is to translate exit angles  $\beta_n$  into corresponding distances from the center of the detector along the spectral direction  $X_n$ . This requires first calculating the distance  $d$  between the camera and the central point of dispersion. To find  $d$ , we can make use of the known values  $X_{n_{edge}}$  and  $\beta_{n_{edge}}$ , as shown in equations 5 and 6.

$$\begin{aligned} X_{n_{edge}}(\lambda_c) &= N_{p1} p_D n_{edge} \\ &= (1)(0.0065mm)(374) = \mathbf{2.431mm} \end{aligned} \quad (5)$$

$$\begin{aligned} d &= \frac{X_{n_{edge}}(\lambda_c)}{\tan(\beta_0(\lambda_c) - \beta_{n_{edge}}(\lambda_c))} \\ &= \frac{2.431mm}{\tan(.1062rad)} = \mathbf{40.74mm} \end{aligned} \quad (6)$$

With all these pieces in place, we can finally calculate the distance from the center of the detector  $X$  as a function of  $n$  and  $\lambda$ , as shown in equation 7. This can then be converted to corresponding pixel number  $P_x(n, \lambda)$  using equation 8.

$$\begin{aligned} X_n(\lambda) &= d \tan(\beta_0(\lambda_c) - \beta_n(\lambda)) \\ &= d \tan(\arcsin(G\lambda_c - \sin(\alpha_0)) - \arcsin(G\lambda - \sin(\alpha_n))) \end{aligned} \quad (7)$$

$$P_x(n, \lambda) = \frac{N_{DX}}{2} + \frac{X_n(\lambda)}{p_D} \quad (8)$$

As an example, we can use this equation to calculate the spectrum location for the edgemost slit.

$$\begin{aligned} P_x(n_{edge}, \lambda_{max}) &= \frac{N_{DX}}{2} + \frac{X_{n_{edge}}(\lambda_{max})}{p_D} \\ &= \frac{N_{DX}}{2} + \frac{d \tan(\arcsin(G\lambda_c - \sin(\alpha_0)) - \arcsin(G\lambda_{max} - \sin(\alpha_{n_{edge}})))}{p_D} \\ &= 1000pix + \frac{(40.74mm) \tan(\arcsin(\frac{600}{mm}(.00055mm) - \sin(.1658rad)) - \arcsin(\frac{600}{mm}(.0007mm) - \sin(.2260rad)))}{.0065mm} \\ &\simeq \mathbf{pixel802} \end{aligned}$$

$$\begin{aligned} P_x(n_{edge}, \lambda_{min}) &= \frac{N_{DX}}{2} + \frac{X_{n_{edge}}(\lambda_{min})}{p_D} \\ &= \frac{N_{DX}}{2} + \frac{d \tan(\arcsin(G\lambda_c - \sin(\alpha_0)) - \arcsin(G\lambda_{min} - \sin(\alpha_{n_{edge}})))}{p_D} \\ &= 1000pix + \frac{(40.74mm) \tan(\arcsin(\frac{600}{mm}(.00055mm) - \sin(.1658rad)) - \arcsin(\frac{600}{mm}(.0004mm) - \sin(.2260rad)))}{.0065mm} \\ &\simeq \mathbf{pixel1945} \end{aligned}$$

Therefore, the outermost slit on the DMD will cast a spectrum on the detector spanning pixels 802 to 1945 in the spectral direction